




# Strategic Business Holdings, Inc.

*...Quality Products, Superior Customer Service, & Support*

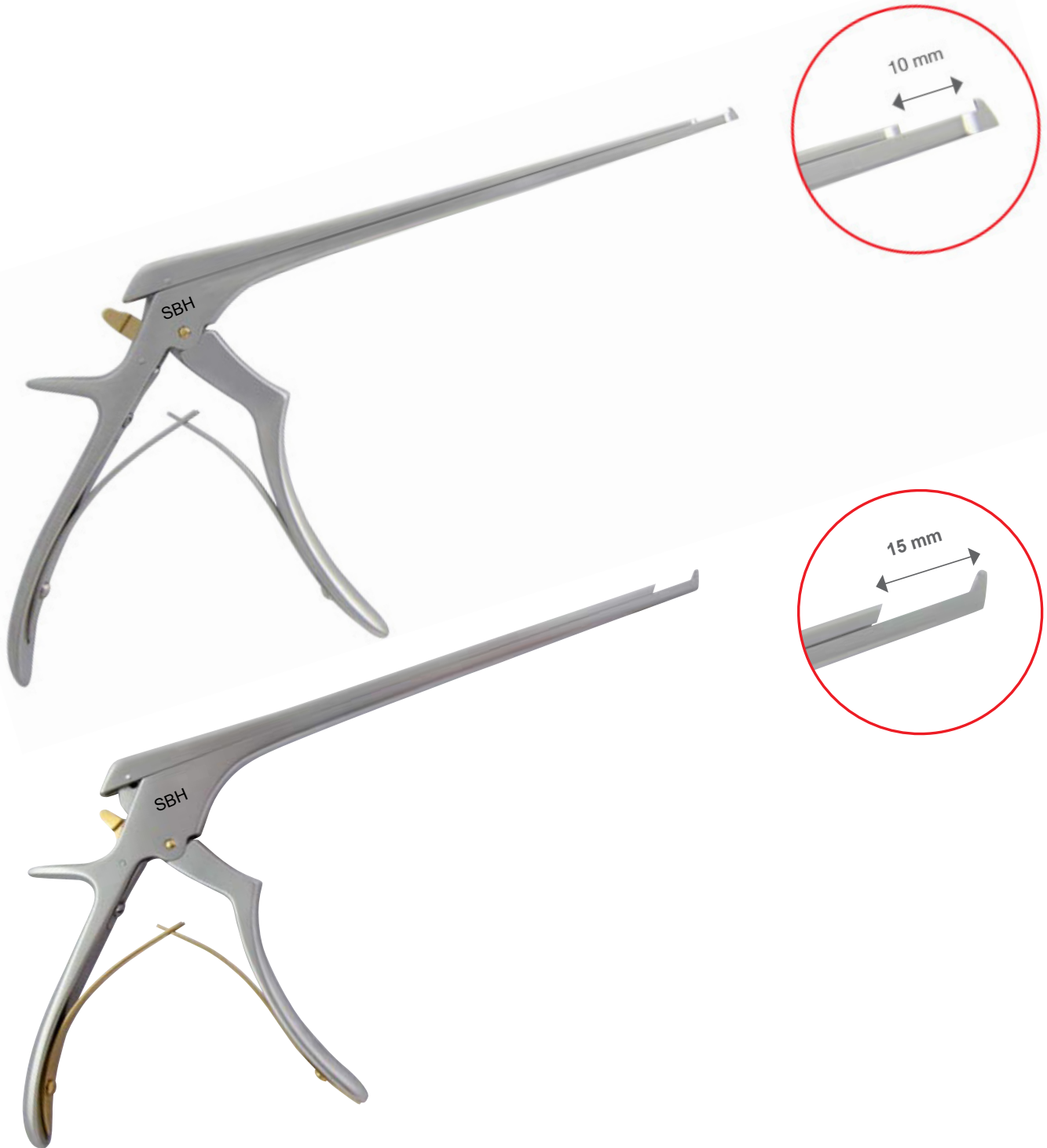

100% made in Germany

## Laminectomy One Piece Detachable Kerrison Rongeur

Copyright © SBH ©. 2020 - All rights reserved.



## Laminectomy Punches Swing-System One Piece Detachable

### Options:

How to order Rongeurs with ejector?

add a "A" to the item number


How to order Rongeurs with thin foot?

add a "T" to the item number

How to order Rongeurs with ceramic coating?


add a "CC" to the item number

Example: You would like to order the S 49-071-01 with ejector and ceramic coating? Order number: S 49-071-01ACC

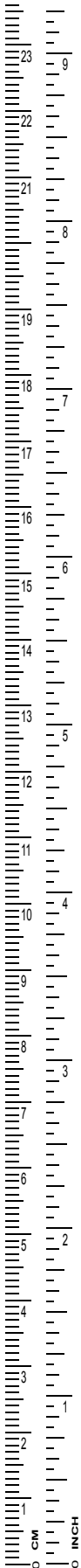



- one piece
- easy to disassemble
- easy to clean
- easy to assemble

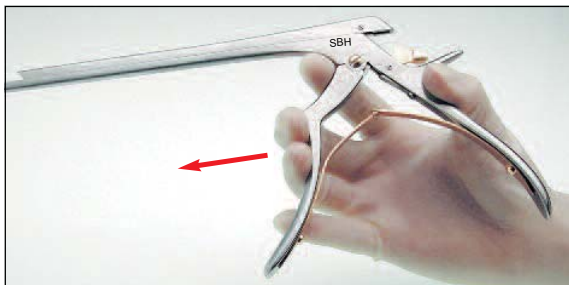


## Disassembly

First close the handle ...

than push the lever downwards ...

now open the handle ...


and push downwards at the end of shaft.

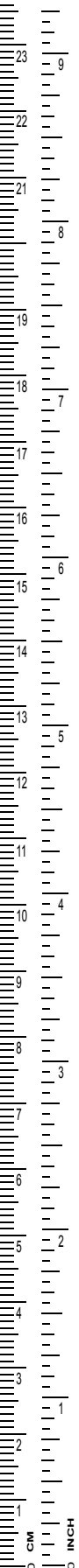



## Laminectomy Punches Swing-System One Piece Detachable

- one piece
- easy to disassemble
- easy to clean
- easy to assemble

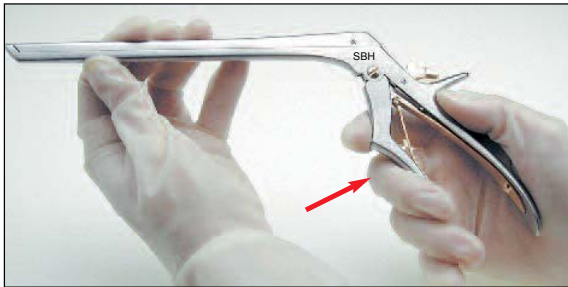


## Assembly

First push shaft parts into one other ...

than close the handle ...

now push the lever upwards ...

and finally open the handle.

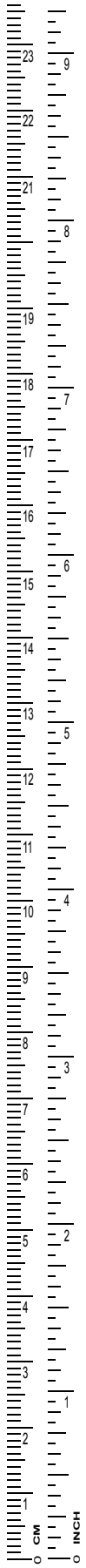


# Laminectomy Punches Swing-System One Piece Detachable



Not Actual Size Shown

## Description


### KERRISON SWING

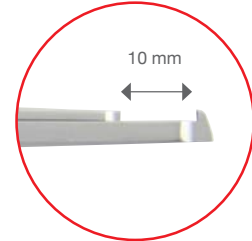

	Width of bite	Length of shaft	Length of shaft	Length of shaft
90°	1 mm	<b>18 cm, 7.2"</b> S 49-051-01 T CC	<b>20 cm, 8"</b> S 49-052-01 T CC	<b>23 cm, 9.2"</b> S 49-053-01 T CC
		without ejector		

	Width of bite	Length of shaft	Length of shaft	Length of shaft
90°	2 mm	<b>18 cm, 7.2"</b> S 49-051-02 T CC	<b>20 cm, 8"</b> S 49-052-02 T CC	<b>23 cm, 9.2"</b> S 49-053-02 T CC
90°	3 mm	S 49-051-03 T CC	S 49-052-03 T CC	S 49-053-03 T CC
90°	4 mm	S 49-051-04 T CC	S 49-052-04 T CC	S 49-053-04 T CC
90°	5 mm	S 49-051-05 T CC	S 49-052-05 T CC	S 49-053-05 T CC
90°	6 mm	S 49-051-06 T CC	S 49-052-06 T CC	S 49-053-06 T CC

with ejector

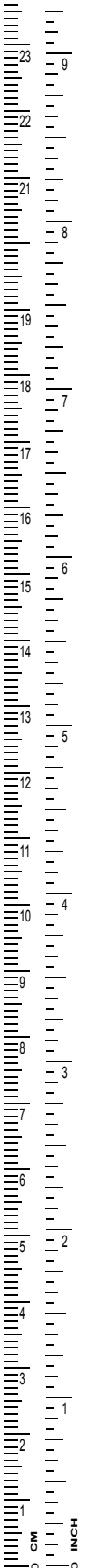




# Laminectomy Punches Swing-System One Piece Detachable

Not Actual Size Shown


**Description**

**KERRISON**

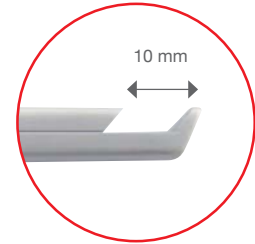

	Width of bite	Length of shaft	Length of shaft	Length of shaft
		<b>18 cm, 7.2"</b>	<b>20 cm, 8"</b>	<b>23 cm, 9.2"</b>
40°	1 mm	S 49-071-01 T CC	S 49-072-01 T CC	S 49-073-01 T CC
40°	2 mm	S 49-071-02 A T CC	S 49-072-02 A T CC	S 49-073-02 A T CC
40°	3 mm	S 49-071-03 A T CC	S 49-072-03 A T CC	S 49-073-03 A T CC
40°	4 mm	S 49-071-04 A T CC	S 49-072-04 A T CC	S 49-073-04 A T CC
40°	5 mm	S 49-071-05 A T CC	S 49-072-05 A T CC	S 49-073-05 A T CC
40°	6 mm	S 49-071-06 A T CC	S 49-072-06 A T CC	S 49-073-06 A T CC

without ejector

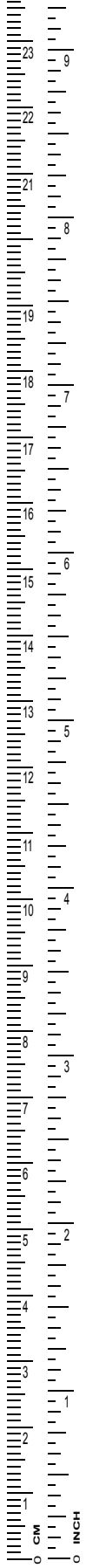


# Laminectomy Punches Swing-System One Piece Detachable



Not Actual Size Shown

**Description**


**KERRISON**

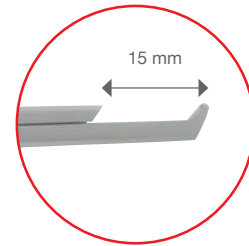

**Width of bite      Length of shaft**

	Width of bite	Length of shaft
		<b>20 cm, 8"</b>
40°	2 mm	S 49-272-02T
40°	3 mm	S 49-272-03T
40°	4 mm	S 49-272-04T
40°	5 mm	S 49-272-05T

without ejector

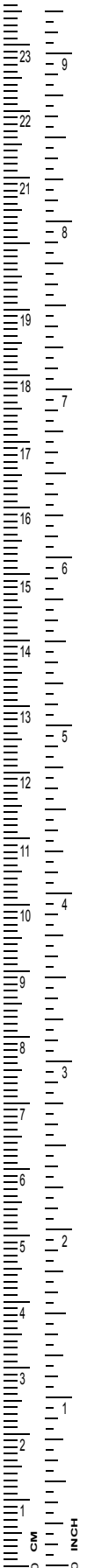



## **Strategic Business Holdings, Inc.**

*...Quality Products, Superior Customer Service, & Support*

Main +1 (800) 303-3050 Fax +1 (949) 600-8874 9620 Research Drive, Irvine - CA 92618 USA [info@sbhsurgical.com](mailto:info@sbhsurgical.com) [www.sbhsurgical.com](http://www.sbhsurgical.com)

Copyright © SBH ®. 2020 - All rights reserved.